

Disconnect Wedge

(en)



DESCRIPTION:

The Mercury Medical® Disconnect Wedge was designed to facilitate disconnections between the 15mm male connectors of endotracheal tubes, tracheostomy tubes or swivel connectors with grasping surfaces (or ears) and the 15mm female connector.


DIRECTIONS FOR USE:

Place the wedge legs (Figure 1) between the male and female connector to be disconnected. While holding the connector closest to the airway, to reduce trauma to the trachea, gently push on the wedge and separate the connections.

NOTE: Never leave a patient with a loosened connection unattended, as an accidental disconnection may result. Repeated disconnection and reconnection may result in an improper and loose fit.

Additional Languages

- (da) Frakobl Wedge
- (nl) Ontkoppel wig
- (fi) Irrota Wedge
- (fr) Déconnecter le coin
- (de) Keil trennen
- (it) Disconnettere il cuneo
- (no) Koble fra kile
- (sr) Дисцоннект Ведре
- (es) Cuña de desconexión
- (sv) Koppla bort Wedge

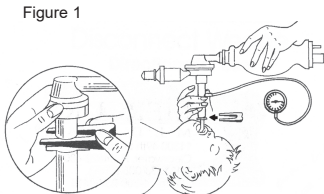

CE
2862



EMERGO EUROPE

Westervoortsedijk 60
6827 AT Arnhem
The Netherlands



Medical
Device



Non
Sterile



Single Patient -
Multiple-Use



Not Made with
Natural Rubber
Latex

Rx

By Prescription
Only