Mercury Medical Superior Solutions
11300 - 49th Street North
Clearwater, Florida 33762-4807 USA
www.mercurymed.com

(en)

FLOW-SAFE II®
DISPOSABLE CPAP SYSTEM

DESCRIPTION:
The Mercury Flow-Safe Continuous Positive Airway Pressure (CPAP) device is a respiratory aid intended for use with a facemask and gas supplying device to elevate pressure in the patient's lungs.

CONNECTIONS:
- Standard oxygen tubing nipple
- Patient connection ISO 5356-1 - 15mm taper female and 22mm taper male

INDICATIONS FOR USE:
The Mercury Flow-Safe CPAP device is intended to provide CPAP to spontaneously breathing patients in the hospital and pre-hospital environment.

WARNINGS:
- Do not allow smoking or use unit near sparking equipment, open flame, oil or other flammable chemicals.
- Any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.

CAUTIONS:
- Federal (USA) law restricts this device to sale by or on the order of a physician.
- Do not clean, soak, rinse or sterilize.
- Reuse of this device may pose a risk of cross-contamination and the device may not perform as intended.
- In the event of undesirable flow rate from oxygen source, simply remove the device and place on supplemental oxygen per protocol.
- Use of the Flow-Safe with non-back pressure compensated flow devices may affect input gas liter flow. Always verify delivered CPAP pressure on manometer.
- This device is assembled with PVC components which include low levels of a phthalate as indicated. Consideration should be used when treating children or pregnant or nursing women.

CONTRAINDICATIONS:
- Respiratory Arrest
- Unconscious
- Cardiogenic Shock
- Pneumothorax
- Facial Anomalies
- Facial Trauma
- Airway Obstruction

DIRECTIONS FOR USE:
1. Connect O_2 tubing nipple to gas source.
2. Secure the face mask snugly to patient's face using head harness.
3. Slowly increase gas flow to 6 or 8 LPM. Check mask fit to patient and device connections for leaks.
4. Adjust the flowmeter until desired pressure is obtained. Flow of 12 - 14 LPM is required to reach CPAP pressure of 8.5 - 10 cm H_2O.
5. Do not exceed 30 LPM.
6. Patient SaO_2 should be monitored using a pulse oximeter.
7. Dispose of according to facility protocol.

MEASURING PRESSURE:
- Pressure relief limits maximum CPAP pressure to 25 cm H_2O @ 25 LPM.
- Do not exceed pressure limit of manometer (25 cm H_2O).
- Manometer accuracy ± 3 cm H_2O up to 15 cm H_2O and ± 5 cm H_2O over 15 cm H_2O.


OPERATING ENVIRONMENTAL LIMITS: -18° C to +50° C

STORAGE ENVIRONMENTAL LIMITS: -40° C to +60° C

MR CONDITIONAL:
Static magnetic field of 3.0 Tesla or less
Spatial Gradient of 3000 gauss/cm or less

CE 2862

EC REP EMERGO EUROPE
Westervoortsedijk 60
6827 AT Arnhem
The Netherlands

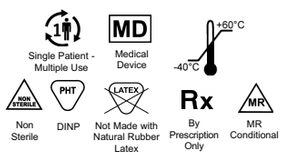® Registered U.S. Trademark of Mercury Enterprises, Inc.
Patent # US 8,522,618 B1; US 9,370,635 B2; US 11,383,054 B2;
US 12,064,551 B2; Other Patents Pending

Mercury Medical Superior Solutions
11300 - 49th Street North
Clearwater, Florida 33762-4807 USA
www.mercurymed.com

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© Εμπορικό σήμα κατατεθέν στις Η.Π.Α. της εταιρείας Mercury Enterprises, Inc.
® U SAD registrovani zaštitni znak Mercury Enterprises, Inc.