Mercury Medical
Superior Solutions
11300 - 49th Street North
Clearwater, Florida 33762-4807 USA
www.mercurymed.com

(en)

FLOW-SAFE II EZ®
DISPOSABLE CPAP SYSTEM

DESCRIPTION:
The Mercury Flow-Safe II EZ Continuous Positive Airway Pressure (CPAP) device is a respiratory aid intended for use with a face mask, nebulizer and gas supplying device to elevate pressure in the patient's lungs while delivering aerosolized medication.

CONNECTIONS:
- Standard oxygen tubing nipple
- Patient connection ISO 5356-1 - 15mm taper female and 22mm taper male
- Nebulizer connection 15mm taper female and 22mm taper male

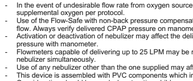
INDICATIONS FOR USE:
The Mercury Flow-Safe II EZ CPAP device is intended to provide CPAP to spontaneously breathing patients in the hospital and pre-hospital environment.

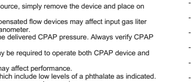
WARNINGS:
- Do not allow smoking or use unit near sparking equipment, open flame, oil or other flammable chemicals.
- Any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.

CAUTIONS:
- Federal (USA) law restricts this device to sale by or on the order of a physician.
- Do not clean, soak, rinse or sterilize.
- Reuse of this device may pose a risk of cross-contamination and the device may not perform as intended.
- In the event of undesirable flow rate from oxygen source, simply remove the device and place on supplemental oxygen per protocol.
- Use of the Flow-Safe with non-back pressure compensated flow devices may affect input gas liter flow. Always verify delivered CPAP pressure on manometer.
- Activation or deactivation of nebulizer may affect the delivered CPAP pressure. Always verify CPAP pressure with manometer.
- Flowmeters capable of delivering up to 25 LPM may be required to operate both CPAP device and nebulizer simultaneously.
- Use of any nebulizer other than the one supplied may affect performance.
- This device is assembled with PVC components which include low levels of a phthalate as indicated. Consideration should be used when treating children or pregnant or nursing women.

CONTRAINDICATIONS:
- Respiratory Arrest
- Unconscious
- Cardiogenic Shock
- Pneumothorax
- Facial Anomalies
- Facial Trauma
- Airway Obstruction

DIRECTIONS FOR USE:
1. Connect O_2 tubing nipple to gas source.
2. Place the face mask securely to the patient's face using head harness.
3. With nebulizer in the OFF position **(Figure 1)**, slowly increase gas flow to 6 or 8 LPM. Check face mask fit to patient and device connections for leaks.
4. Adjust the flowmeter until desired pressure is obtained. Flow of 12 - 14 LPM is required to reach CPAP pressure of 8.5 - 10 mm H_2O.
5. Do not exceed 30 LPM.
6. Patient SaO_2 should be monitored using a pulse oximeter.
7. To activate nebulizer rotate knob to ON position **(Figure 2)**.
8. If necessary, readjust flowmeter to obtain desired CPAP pressure. Up to 25 LPM may be required.
9. Dispose of according to facility protocol.

MEASURING PRESSURE:
- Pressure relief limits maximum CPAP pressure to 25 cm H_2O @ 25 LPM.
- Do not exceed pressure limit of manometer (25 cm H_2O).
- Manometer accuracy ± 3 cm H_2O up to 15 cm H_2O and ± 5 cm H_2O over 15 cm H_2O.

SPECIFICATIONS:
Sample Guidelines for Preparing RX Dosing

FLOWMETER SETTING L/min	14 – 15				23 – 24			
CPAP PRESSURE cm H_2O	4 – 5				9 – 10			
FLOW THROUGH EZFLOW MAX	6 L/min				10 L/min			
OUTPUT	12 mL/hour				16 mL/hour			
RX (mg/hour)	5	10	15	20	5	10	15	20
Treatment Duration (hours)	1 2	1 2	1 2	1 2	1 2	1 2	1 1.5	1 1.5
Medication @ 5mg/mL (mL)	1 2	2 4	3 6	4 8	1 1.5	2 3	3 4.5	4 6
Saline (mL)*	11 22	10 20	9 18	8 16	15 22	14 21	13 20	12 18

** Rounded to the nearest mL*

EXAMPLE:
To deliver 10 mg/hour of 5 mg/mL medication at 10 L/min; *1.5 Hour Duration:* **Place 3 mL of medication + 21 mL of Saline** into nebulizer reservoir and run for 1.5 hours.

CAUTION: *The listed output is nominal value only, actual output may hervivier kun vary depending on device, accuracy of flowmeter, and the existence of tubing connection leaks. Follow up output checks are recommended and flowmeter adjustments may be necessary.*

OPERATING ENVIRONMENTAL LIMITS: -18° C to +50° C
STORAGE ENVIRONMENTAL LIMITS: -40° C to +60° C

MR CONDITIONAL:
Static magnetic field of 3.0 Tesla or less
Spatial Gradient of 3000 gauss/cm or less

CE 2862
EC REP **EMERGO EUROPE**
Westervoortsedijk 60
6827 AT Arnhem
The Netherlands

® Registered U.S. Trademark of Mercury Enterprises, Inc.
Patent # US 8,522,618 B1, US 6,338,443 B1, US 9,370,635 B2, US 9,511,202,
US 11,383,054 B2; US 12,064,551 B2; Other Patents Pending

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© Gedeponeerd Amerikaans handelsmerk van Mercury Enterprises, Inc.
® Marca de Mercury Enterprises, Inc. registrada en los EE.UU.