

AmbuMan

AmbuMan is the ideal training manikin for BLS CPR. It is specifically designed for the basic training of medical and nursing staff as well as lay people.

AmbuMan is designed to provide optimal training of chest compressions and demonstrate realistic respiration. The user friendly Ambu hygienic system prevents cross-infection and provides realistic airway resistance and exhalation through mouth and nose ventilation.

Monitoring Instrument

The indicator shows the instructor the volume of supplied air (in L) and the chest compression depth (in mm). Furthermore, a warning symbol is shown when a trainee is performing compressions with incorrect hand placement and when there is gastric inflation.

Carotid Pulse

The carotid pulse can be felt during chest compression and can be simulated by a manual pulse.

AmbuMan® Instrument



*Arms and legs are optionally available

KEY BENEFITS

AmbuMan Instrument

Optimal BLS manikin

Robust design

Realistic training of chest compressions

Adjustable chest stiffness

Ambu hygienic system for mouth-to-mouth and mouth-to-nose ventilation

Airways open only with correct hyperextension


Manual feedback display (ventilation volume, compression depth, hand position, and stomach inflation)

The carotid pulse can be felt during chest compressions and can also be simulated manually

Ambu

SPECIFICATIONS

Connectivity	WLAN or LAN Ambu Wireless Technology built on the latest generation of Intel WLAN Modules
Ambu Manikin Management Module	Operating system independent and browser-based control, programmable scenarios, individual documentation for up to 4 participants, parallel display and control of up to 6 manikins, parallel multiple access of various devices for visualization of the performance
Power supply	AC adapter, Battery pack or PoE (optional, not included in packaging)
Manual feedback display	Extendable Ventilation volume, compression depth, hand position, stomach inflation

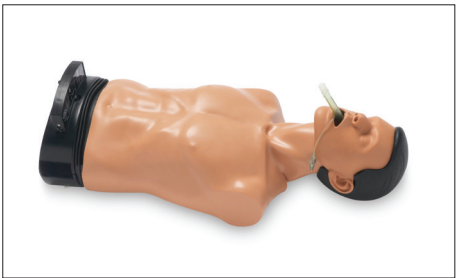

*Arms and legs are optionally available for the AmbuMan Instrument

TECHNICAL DATA

Weight torso with transport bag	12 kg
Length torso	80 cm
Length full body	165 cm



Monitoring instrument

ORDER INFORMATION

Item No.	Description
A 234 407 000	AmbuMan Wireless Torso, Integrated Ambu Wireless Technology, Ambu Manikin Management Module, 5 Ambu face pieces, 100 disposable head bags, Transport bag/practice mat
A 234 404 000	AmbuMan Instrument Torso, 5 Ambu face pieces, 100 disposable head bags, Transport bag / practice mat
234 000 077	Battery pack
234 000 706	Upgrade kit full body
234 000 702	Ambu head bags, disposable (set of 100)
234 000 703	Ambu face pieces (set of 5)
234 000 725	Ambu Guedel face pieces (set of 5)



Transport bag (practice mat incl.)

Questions? Contact your Ambu representative at **800-262-8462** or visit **www.ambuUSA.com/**



Ambu Inc.
6721 Columbia Gateway Drive
Suite 200
Columbia, MD 21046
Tel. 800 262 8462
Fax 800 262 8673
ambuUSA.com

