

AIR-Q® sP3G

SUPRAGLOTTIC AIRWAY



AIR-Q®sp3G

INTUBATING LARYNGEAL AIRWAYS (ILA)




The advanced, all-silicone Air-Q®sp3G helps deliver safety and confidence in airway management, no matter what airway challenges you face. The unique ET tube ramp and epiglottis elevator aid in fast, easy and safe intubations.

SPECIFICALLY DESIGNED FOR EMS IN MIND

Air-Q®sp3G - Features

- Uses positive pressure ventilation to self-pressurize the cuff
- Pressures cycle between peak airway pressure and PEEP levels
- Has an extra wide gastric inlet access and two gastric channels that accommodate NG tubes up to 18 Fr

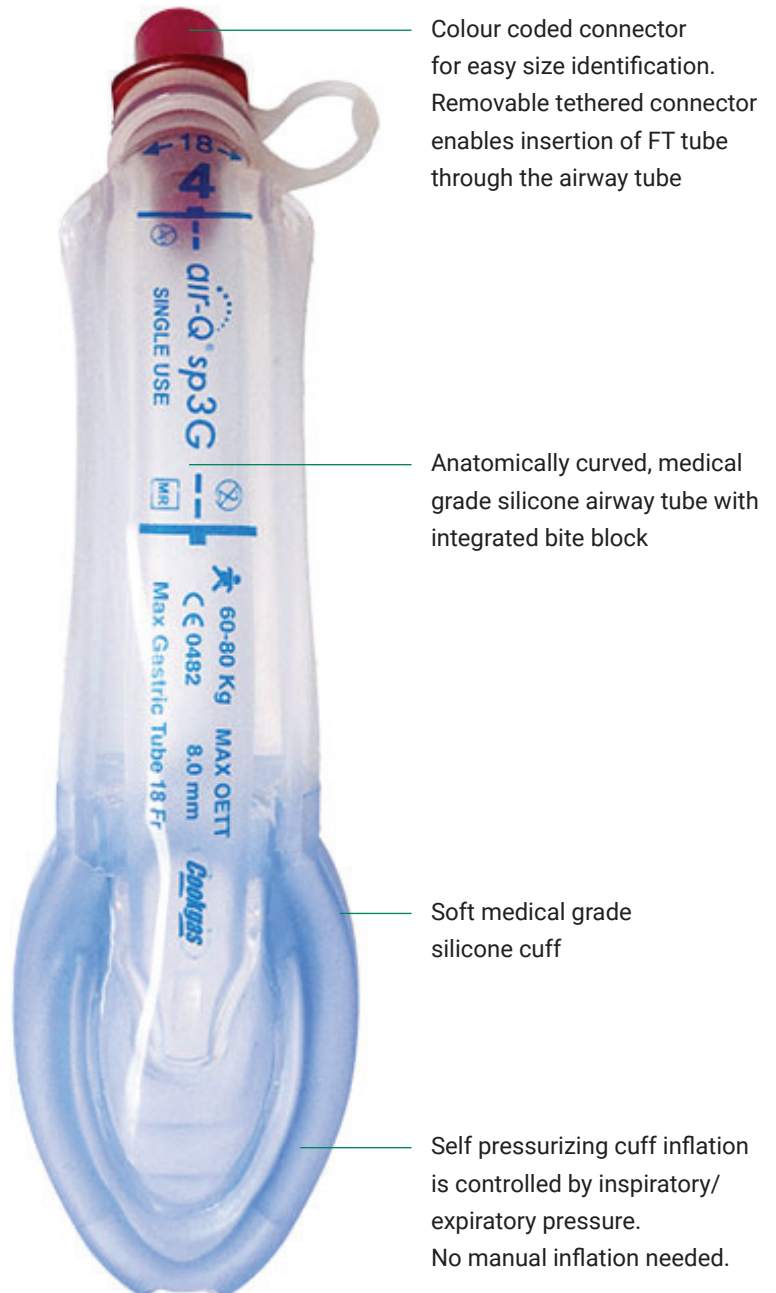

AIR-Q®sp3G

PERFECT PLACEMENT MADE EASY




Air-Q®sp3G


SUPRAGLOTTIC AIRWAY PROVIDES

- Anatomically correct placement of the masked airway
- Lifting of the epiglottis to maximize access
- Guide ramp to assist with ETT insertion and placement

Air-Q®sp3G PROVIDES BETTER VISUAL CLARITY

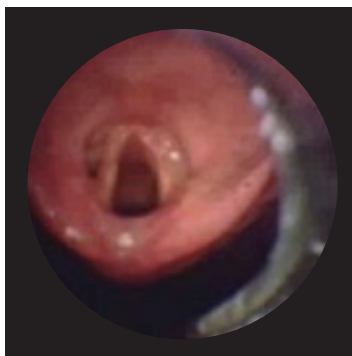


Air-Q®sp3G Series
with Epiglottis Elevator

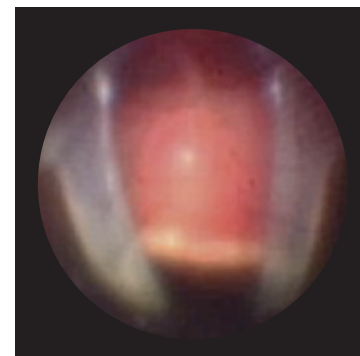

Competitive ILA

EPIGLOTTIS ELEVATOR

Our unique supraglottic airway includes an epiglottis elevator and raised heel to optimize airway placement and seal.

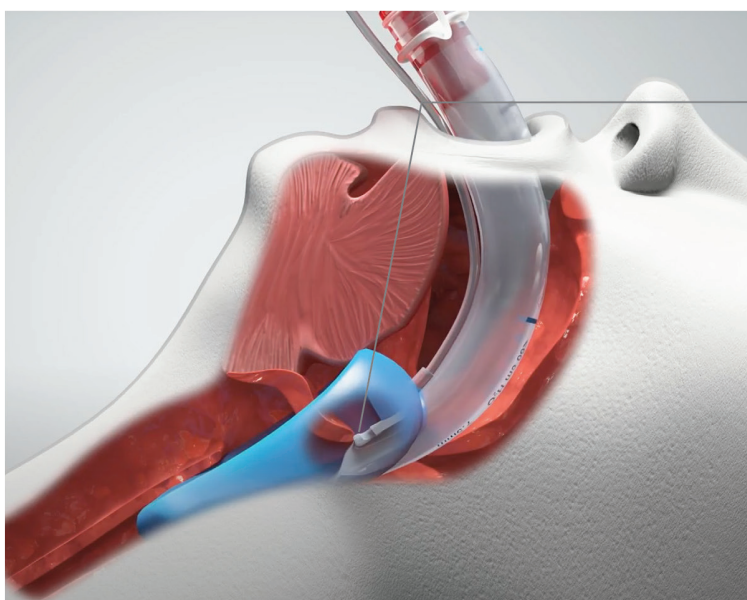

Built-up mask heel for improved seal

Epiglottis Elevator

ET tube ramp helps facilitate intubation and directs the endotracheal tube toward the laryngeal inlet



AIR-Q®sp3G

PERFECT PLACEMENT MADE EASY



AIR-Q®sp3G INTUBATING LARYNGEAL AIRWAY WITH 15MM CONNECTOR

ITEM	DESCRIPTION	SIZE	ADAPTER COLOR	IDEAL BODY WEIGHT	MAXIMUM OEETT	MOUTH OPENING	LENGTH OF TUBE	MAXIMUM OG TUBE	ET TUBE CONNECTOR
2114-60005	Air-Q®sp3Q	0	Light Purple	<2.0kb	3.0 mm	5 mm	6.0 cm	5 Fr	15 mm
2114-60055	Air-Q®sp3Q	0.5	Pink	2 – 4 kg	4.0 mm	8 mm	7.0 cm	6 Fr	15 mm
2114-60105	Air-Q®sp3Q	1.0	Blue	4 – 7 kg	4.5 mm	11 mm	9.0 cm	8 Fr	15 mm
2114-60155	Air-Q®sp3Q	1.5	Green	7 – 17 kg	5.0 mm	14 mm	11.0 cm	8 Fr	15 mm
2114-60205	Air-Q®sp3Q	2	Orange	17 – 30 kg	5.5 mm	17 mm	14.0 cm	10 Fr	15 mm
2114-60305	Air-Q®sp3Q	3	Yellow	30 – 60 kg	7.0 mm	20 mm	16.0 cm	14 Fr	15 mm
2114-60405	Air-Q®sp3Q	4	Red	60 – 80 kg	8.0 mm	23 mm	18.0 cm	16 Fr	15 mm
2114-60505	Air-Q®sp3Q	5	Purple	> 80 kg	9.0 mm	25 mm	20.0 cm	18 Fr	15 mm