

CuraSafe™ Comfort Push Button Safety IV Catheter

VenaPUR® (polyurethane) Material designed to remain kink-free


product specifications

- Push button style needle retraction
- Ground stainless-steel needle
- Straight hub
- Disposable
- Non-pyrogenic
- Not made with natural rubber latex
- Not made with DEHP
- ISO 10555-5 compliant color-coding
- Needlestick Safety and Prevention Act (NSPA) of 2000 compliant

Triple-Facet, Back Cut, Beveled Needle
Designed to provide ease of insertion

Gently Tapered Catheter Tip
Provides smooth transition from needle to catheter during insertion


Push Button Style Needle Retraction Safety Mechanism
Retracts needle fully and safely into catheter body

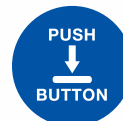

Flashback Chamber
Allows for clear visual confirmation of successful venipuncture

Push-Off Tab
Assists with one-handed catheter advancement

Encapsulating Needle Shield
Encases the needle when safety mechanism is activated, promoting needlestick safety and reducing exposure to bloodborne pathogens

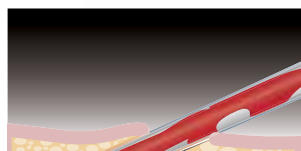










Featured in sizes 20G - 24G

RapidFlash® Notched Needle
Allows for quick, clear visual confirmation of successful venipuncture

Ridged Catheter Hub
Offers superior grip and stability during connection of IV Administration Sets

Contoured Catheter Body
Ergonomically designed for one-handed technique

ordering information

Color	Item #	Gauge and Length	Catheter ID (mm)	Catheter OD (mm)	Gravity Flow Rate (mL/min)	Qty
	1612-84310	14G x 1.75"	1.6mm	2.1mm	297mL/min	50/bx, 200/case
	1612-84320	16G x 1.77"	1.38mm	1.8mm	220mL/min	50/bx, 200/case
	1612-84330	18G x 1.88"	0.95mm	1.3mm	94mL/min	50/bx, 200/case
	1612-84340	18G x 1.16"	0.95mm	1.3mm	94mL/min	50/bx, 200/case
	1612-84360	20G x 1.88"	0.78mm	1.1mm	61mL/min	50/bx, 200/case
	1612-84350	20G x 1.16"	0.78mm	1.1mm	61mL/min	50/bx, 200/case
	1612-843203	20G x 1"	0.78mm	1.1mm	61mL/min	50/bx, 200/case
	1612-84370	22G x 1"	0.62mm	0.9mm	36mL/min	50/bx, 200/case
	1612-84380	24G x 0.75"	0.49mm	0.7mm	20mL/min	50/bx, 200/case

